



# Product Specifications


Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

## Oligo Modifications

For research use only. Not for use in diagnostic procedures for clinical purposes.

### Cy5.5

Category	Fluorescent Dyes
Modification Code	Cy5.5-5
Reference Catalog Number	26-6460
5 Prime	Y
3 Prime	N
Internal	N
Molecular Weight(mw)	633.74

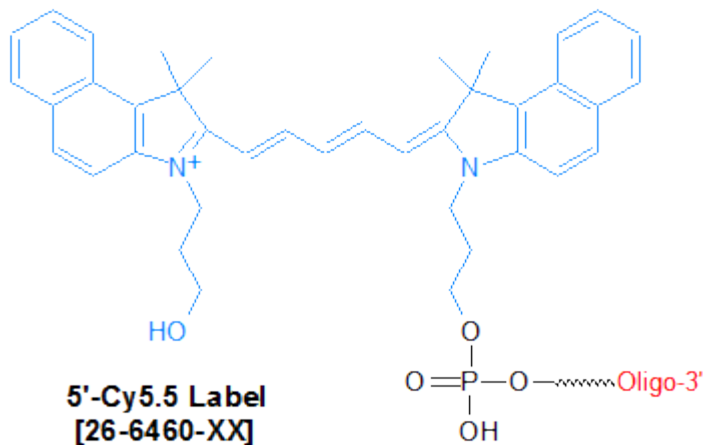

Cyanine 5.5 (Cy5.5) is a far-red/near infra red (NIR) fluorescent dye that belongs to the Cyanine family of synthetic polymethine dyes. Cy5.5 is reactive, water-soluble, and has an absorbance maximum of 675 nm and an emission maximum of 694 nm. It is available as a phosphoramidite, and is used to fluorescently label oligonucleotides at either the 5' or 3' end, or internally. Cy5.5 plays a particularly important role in real-time PCR applications, being used as a reporter moiety in TaqMan probes (1), Scorpion primers (2) and Molecular Beacons (3). For such probes, Cy5.5 is most commonly paired with the dark quencher BHQ-3, as the two have excellent spectral overlap.

Cy5.5 can also be used to label DNA oligos for use as hybridization probes in other applications, such as Fluorescent In-Situ Hybridization (FISH).

Near Infrared Fluorophore Spectral Data & Quencher Selection Guide

Fluorophore Name

Excitation Max, nm +/-10

Emission Max, nm +/-10

Extinction Coefficient\*

Color\*\*

Quencher

**Cy5 650 665 250,000**

**IRDye 650 NHS 650 665 230,000**

asp?modid=516">AZ647 NHS 655 680 191,800

**Cy5.5 684 710 198,000**

**IRDye 700 NHS 684 710 288,000**

**Cy7 NHS 740 773 199,000**

**IRDye 750 NHS 756 776 260,000**

**cy7.5 NHS 788 808 223,000**

**IRDye 800 NHS 795 819 240,000**

**\* Extinction coefficient at  $\lambda$  (max) in  $\text{cm}^{-1}\text{M}^{-1}$ . \*\* Typical emission color seen through the eyepiece of a conventional fluorescence microscope with appropriate filters. Near-IR region. Human vision is insensitive to light beyond ~650 nm; it is not possible to view near-IR fluorescent dyes.**

**[Click here for a list of fluorophores.](#)**

[Click here for list of quenchers.](#)

#### References

1. Livak, K.J., Flood, S.J.A., Marmaro, J., Giusti, W., Deetz, K. Oligonucleotides with fluorescent dyes at opposite ends provide a quenched probe system useful for detecting PCR product and nucleic acid hybridization. *PCR Methods Appl.* (1995), 4: 1-6.
2. Thelwell, N., Millington, S., Solinas, A., Booth, J., Brown, T. Mode of action and application of Scorpion primers to mutation detection. *Nucleic Acids Res.* (2000), 28: 3752-3761.
3. Tyagi, S., Kramer, F.R. Molecular beacons: probes that fluoresce upon hybridization. *Nat. Biotechnol.* (1996), 14: 303-308.
4. Stoecker, K., Dorninger, C., Daims, H., Wagner, M. Double Labeling of Oligonucleotide Probes for Fluorescence In Situ Hybridization (DOPE-FISH) Improves Signal Intensity and Increases rRNA Accessibility. *Appl. Environ. Microb.* (2010), 76: 922-926.