



## Product Specifications


Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

## Oligo Modifications

For research use only. Not for use in diagnostic procedures for clinical purposes.

### 5'-BBQ-650

Category	Quenchers
Modification Code	5-BBQ-650
Reference Catalog Number	26-6696
5 Prime	Y
3 Prime	Y
Internal	Y
Molecular Weight(mw)	667.63

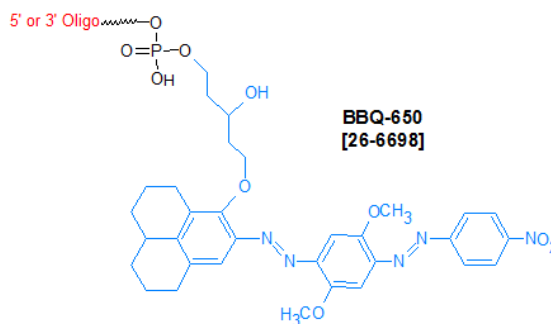

BlackBerry Quencher 650 (BBQ650) is classified as a dark quencher (a non-fluorescent chromophore). Dark quenchers are extensively used as the 3'-quencher moiety in a variety of Fluorescence Resonance Energy Transfer (FRET) DNA detection probes in which the fluorophore has a long wavelength (yellow to far red) emission maximum (e.g. Cy3, ROX, Cy5, Cy 5.5). Dark quenchers can serve in this role because they have long wavelength absorbance maxima. Dark quenchers are primarily used in nucleic acid assays, but also find a place in nucleic acid structural studies (1). Examples include TaqMan probes (2), Scorpion primers (3), and Molecular Beacons (4).

BBQ650 has an absorbance maximum of 650 nm, and an effective absorbance range of 550-750 nm (yellow to far red). It is chemically resistant to both oligonucleotide synthesis reagents (iodine, TCA) or deblocking solutions (ammonia, AMA).

**Consequently, for synthesis of longer oligos (> 50 bases), BBQ650 is the preferred quencher over BHQ-2 or BHQ-3, as the latter are chemically less stable, and degrade when exposed to oligo synthesis and deprotection conditions for long periods of time, such as when synthesizing or processing longer oligos.**

A list of specific fluorescent dyes compatible with BBQ650 is found at this link.

[Click here for complete list of quenchers](#)

The emission spectra of this set of dyes sufficiently overlaps the absorbance spectrum of BBQ650 to allow the latter to quench the fluorescence of the former with a high degree of efficiency.

The advantages of using BBQ650 as a dark quencher in a FRET probe are (a) low background fluorescence (and thus better signal-to-noise ratio), (b) higher dynamic range, (c) amenability to multiplex assays (due to a dark quencher having no secondary fluorescence), and (d) ease of synthesis of FRET probes (due to BBQ650 being resistant to degradation during the oligo deprotection step). <

[Click here for list of quenchers.](#)

[Click here for a list of fluorophores.](#)

Quencher Spectral Data

Quencher

Absorption Max, nm

Quenching Range, nm Dabcyl 453 380-530 BHQ-0 495 430-520 BHQ1 534 480-580 BHQ2 579 550-650 BHQ3 672 620-730  
BBQ-650 650 550-750

[genelink.com/oligo\\_modifications\\_reference/OMR\\_mod\\_category\\_intro.asp?mod\\_sp\\_cat\\_id=15](http://genelink.com/oligo_modifications_reference/OMR_mod_category_intro.asp?mod_sp_cat_id=15)>Click here for complete list of quenchers and details \*\*Black Hole Quencher License Agreement

Black Hole Quencher License Agreement. "Black Hole Quencher<sup>®</sup>, BHQ<sup>®</sup>, CAL Fluor<sup>®</sup> and Quasar<sup>®</sup> are registered trademarks of Biosearch Technologies, Inc., Petaluma, California. The BHQ, CAL Fluor and Quasar dye technologies are protected by U.S. and world-wide patents either issued or in application. Compounds incorporating these dyes are made and sold under agreement with Biosearch Technologies, Inc. for end-user's non-commercial research and development use only. Their use in commercial applications is encouraged but requires a separate Commercial Use License granted by Biosearch Technologies, Inc."

#### References

1. Didenko, V.V. DNA Probes Using Fluorescence Resonance Energy Transfer (FRET): Designs and Applications. *Biotechniques* (2001), **31**: 1106-1121.
2. Livak, K.J., Flood, S.J.A., Marmaro, J., Giusti, W., Deetz, K. Oligonucleotides with fluorescent dyes at opposite ends provide a quenched probe system useful for detecting PCR product and nucleic acid hybridization. *PCR Methods Appl.* (1995), **4**: 1-6.
3. Thelwell, N., Millington, S., Solinas, A., Booth, J., Brown, T. Mode of action and application of Scorpion primers to mutation detection. *Nucleic Acids Res.* (2000), **28**: 3752-3761.
4. Tyagi, S., Kramer, F.R. Molecular beacons: probes that fluoresce upon hybridization. *Nat. Biotechnol.* (1996), **14**: 303-308.